





INTRAORAL IMAGING PLATE SCANNER SO COMPACT, YET SO EFFECTIVE...



OWANDY-CR²

THE MOST COMPACT AND EFFECTIVE INTRAORAL IMAGING PLATE SCANNER... ONCE YOU'VE TRIED IT, YOU WON'T WANT TO USE ANYTHING ELSE!

Optimised sensitivity

The **Owandy-CR**² features a special cuttingedge concept reading head that optimises the way in which plates are scanned. Using this tool, you can produce **accurate** images at **low doses** without the risk of overexposure. The **Owandy-CR**² incorporates specially designed image filters and processing to enable you to produce quality images with excellent contrasts.



DEVICE IS SUITABLE TO ALL YOUR CLINICAL PRACTICE



Endodontics



Prosthetic and implant surgery



Periodontics



Owandy

Caries diagnosis





Diminutively-proportioned

The **Owandy-CR**² is the smallest system on the market! It will easily find its place in all dental practices: its compact, stylish design will add a contemporary and hi-tech touch to your treatment room!



100% automatic & safe

The **Owandy-CR**² is fully automatic! Getting to grips with the system is both quick and easy.

The **Owandy-CR**² is highly intuitive: imaging plate detection, their size and their insertion direction, reading and erasing are automatic: **no complicated user interface**. The unit comprises an **automatic wake up and standby mode**.

For 100% reliable hygiene, the plates are used with cardboard protection and a singleuse hygienic envelope.







OWANDY-CR²

SO COMPACT, YET SO EFFECTIVE YOU'LL NEVER WANT TO BE WITHOUT IT!



Suitable for all types of clinical practice

The **Owandy-CR**² can be used for endodontics, prosthetic dentistry, implant surgery and periodontics, as well as in caries diagnosis.

Different plate sizes are available to meet your requirements (bitewing, periapical, posterior images, etc.) and your patients:

- Size 0: Child
- Size 1: Incisor block
- Size 2: Adult
- Size 3: Bitewing
- Size 4: (2 x size 3): Occlusal block




Ethernet or USB

Whether you opt for the USB or Ethernet version, the **Owandy-CR²** makes it easier to use and install than any other model. Once the imaging plate has been scanned, the image taken can instantly be viewed on the computer screen!

Whether USB or Ethernet, the **Owandy-CR²** has been developed for a single or multi-user use. Sharing data becomes easier!





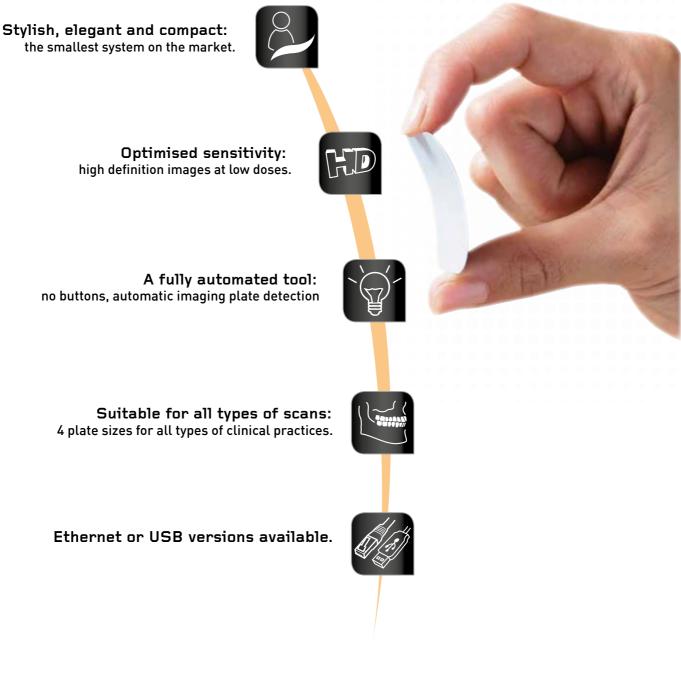





Owandy Radiology is an internationally-recognised manufacturer of dental radiology systems and developer of imaging software. Its software is currently used by almost half of the dental practitioners in France and its range of digital imaging solutions is made available to dentists in more than 50 countries worldwide. In addition to its products and services, **Owandy Radiology** strives to contribute to the improvement of oral health and patients' well-being. As a committed and innovative company, **Owandy Radiology** has been a partner of choice for dentists for more than 20 years.

TECHNICAL FEATURES

Technology	Re-usable imaging plate PSP: Photostimulable Phosphor plate
Resolution	34 pixel/mm – 16 Bits
Pixel size	30 µm
Scan time	4 - 9 s
Connection	USB 2.0 & USB 3.0
	Ethernet 10/100/1000 Mbps
Dimensions	286 x 151 x 104 mm
Weight	4.25 kg
External power supply	100-240 V, 50/60Hz, 60W
Plate size selection	Automatic detection
System standby mode	Automatic function
Plate sizes	Size 0:31 x 22mm, 1024 x 726 pixels
	Size1:40 x 24mm, 1321 x 792 pixels
	Size 2 : 41 x 31mm, 1354 x 1024 pixels
	Size 3:54 x 27mm, 1783 x 891 pixels
Operating system	Windows 7, 8, 10 (32-64 Bits)

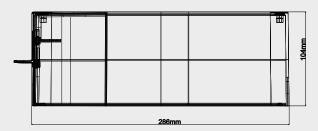

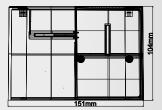

Medical device class IIa/CE0051. Read carefully the instructions before use.

Dealer's stamp	ealer's stamp)
----------------	---------------	---

Owandy Radiology

2, rue des Vieilles Vignes 77183 Croissy-Beaubourg FRANCE Tel.: +33 (0)1 64 11 18 18 info@owandy.com www.owandy.com



DIGITALSENSORS

IMAGING PLATE SCANNER 2D - 3D PANORAMICS INTRAORAL RADIOLOGY CAMERAS SOFTWARE