




**MEDVAL**  
healthcare

**ips**

*Innovation for better life*


# A new era of healthcare

Dermatology Portfolio

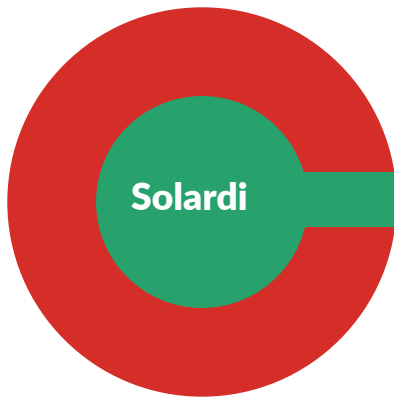

<b>Laser Source</b>	
<b>Laser Type</b>	Diode Laser
<b>Wave lenght</b>	808 nm
<b>Energy</b>	1-60 J/cm adjustable each 1J/cm
<b>Function</b>	
<b>Operating Mode</b>	Arm, Bread, Breast, Armpit, Groin, Leg
<b>Spot Speed</b>	1-15 Hz
<b>Puse Width</b>	7.5- 1200 ms (600W) 10-2025 ms (300W)
<b>Beam Shape</b>	11 x15mm
<b>Light Conduction</b>	Crystal Coupling
<b>Cooling Tip Temperature</b>	Up to - 5°C
<b>General Features</b>	
<b>Power Supply</b>	220-240 V~, 50/60Hz
<b>Cooling System</b>	Water Cooling equipped with Chiller
<b>Dimension</b>	500(W) x600(D)x1470(H)mm
<b>Weight</b>	75 Kgs
<b>Display</b>	Touch Screen Color LCD

## How does it work

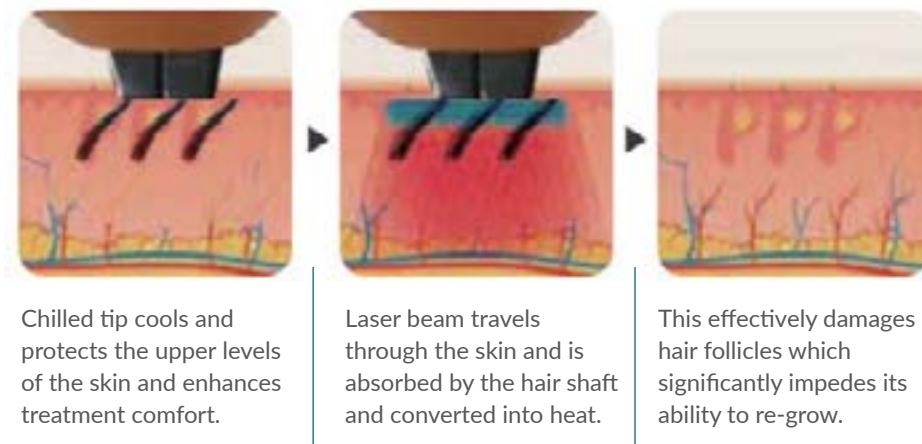

Chilled tip cools and protects the upper levels of the skin and enhances treatment comfort.

Laser beam travels through the skin and is absorbed by the hair shaft and converted into heat.

This effectively damages hair follicles which significantly impedes its ability to re-grow.

### Features:

- Fits all skin types and hair types.
- Various pulse modes: Short, long and combi.
- Pre-set parameters.
- No consumables.
- Crystal contact cooling up to -5 degrees C.



Before

After



Before

After





**MEDVAL**  
healthcare

**ips**

*Innovation for better life*



Villa 8 block 5, Al-Mokattam club  
street, El-Togarieen city, Al-Mokattam,  
Cairo, Egypt.



+202 25086315



+202 25086310

[www.medval-emea.com](http://www.medval-emea.com)

